

MAXIMUM SECURITY



FORMALDEHYDE DISPENSER

Model DSF



AnatHomic Solutions
PATHOLOGY LABORATORY

Why to choose our Formaldehyde Dispenser?

- ✓ Because it offers maximum safety and comfort to the user: it has a powerful rear and upper ventilation, and ventilation in the formaldehyde storage cubicle.
- ✓ Because it is possible to configure it according to the needs of the user
- ✓ Because its glossy finish facilitates optimal sanitation
- ✓ Because the walls are covered with insulating material to reduce noise as much as possible
- ✓ Because it is the only dispenser on the market that includes the OZONIC Wash System® to disinfect wastewater
- ✓ Because it has a filter counter that allows to maximize its use
- ✓ Because it is the only dispenser on the market that uses biodegradable filters and very easy to replace
- ✓ For its convenient replacement price of filters manufactured by Anathomic Solutions



GENERAL CHARACTERISTICS

Stainless Steel construction	AISI – 304
Sheet thicknesses	1,5 - 2 mm
Overall dimensions	1500 x 780 x (h)1990 mm
Connections:	
<ul style="list-style-type: none"> • Single-phase electrical connection: 230V 50Hz + ground • Hot water supply 3/8" • Cold water supply 3/8" • Drain PVC Ø40 mm 	
Electrical power	750 W
Control panel	✓
Safety switch (power sectioning)	✓
Emergency push-button (in case of formaldehyde spills)	✓
Sides with window (safety glass)	Laminar glass 3+3
Front partition with electric guillotine system	Laminar glass 4+4
Levelling feet	✓
Formaldehyde dispenser tap	✓
Lower suction enclosure for storage formaldehyde tanks:	
<ul style="list-style-type: none"> • 1 new formaldehyde deposit 10L • 1 used formaldehyde deposit 25L • Filling and emptying level sensors • Quick connect valves + drip 	✓
Automatic formaldehyde dispenser pedal	✓
Single-lever medical tap (H/C water)	✓
Perforated removable trays	1
Top LED lighting	3
Inlet water pressure	3-8 bar
Functioning temperature	5°C – 45°C

ASPIRATION SYSTEM

Double suction zone: front and top	✓
Centrifugal fan motors	3
Activated carbon filters with potassium permanganate impregnation	3
Filter life control	✓
First prior notice to change filters	800 h
End of life of filters	850 h
Adjustable suction speed 10% up to 100%	Intervals ± 5 %
Maximum suction capacity	1404 m ³ / h
Suction speed in cabin	0.5 - 1.15 m/s
Front screen lowered, distance to work area	125 mm
Front screen up, distance to work area	450 mm

BIODEGRADABLE ACTIVE CARBON FILTERS IMPREGNATED WITH POTASSIUM PERMANGANATE

All our stations are equipped with biodegradable filters made from corn and filled with Active Carbon impregnated with Potassium Permanganate, the most effective way to disinfect exhaled particles into the environment. These filters are manufactured by AnatHomic Solutions.

A unique filter, the most effective, easily replaceable, with the lowest environmental impact. Innovative design and in-house manufacturing from AnatHomic Solutions.




OZONE DISINFECTION SYSTEM

All the AnatHomic Solutions stations are equipped with OZONIC WASH 4.0, which through the power of ozone to eliminate pesticides, fungi, organic matter, pollutants, and viruses. In addition, it is 100% free of toxins, without adding contaminants to the water or the environment.



DIMENSIONS

DSF - 003	1500 mm	780 mm	1990 mm
WASHING SINK	600 mm	450 mm	250 mm
DISPENSING SINK	600 mm	450 mm	150 mm

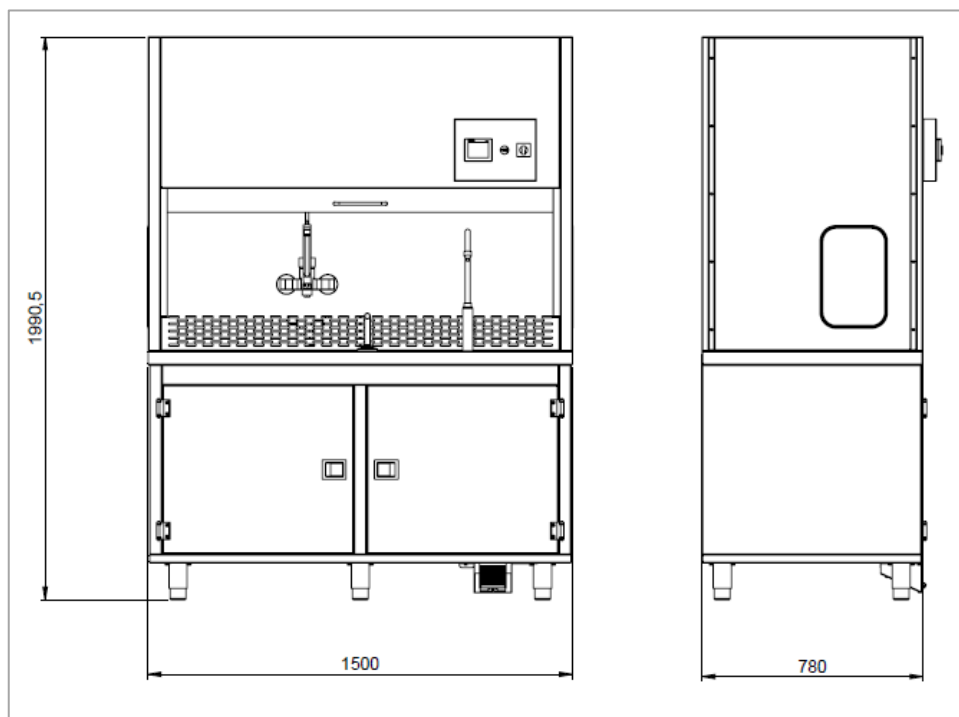

PACKAGING

12-01 Phytosanitary packaging + VCI cover	✓
Wooden box with Anathomic Solutions logo	✓
Number of parcels	1
Dimensions equipment + packaging	
Weight equipment + packaging	



ANATHOMIC SOLUTIONS CERTIFICATIONS



UNE-EN ISO 9001:2015



UNE-EN ISO 13485:2018



UNE-EN ISO 14001:2015



ANATHOMIC SOLUTIONS SECURITY PATHOLOGY LAB CONCEPT

Your Safety is our Motivation for Design and Manufacture Innovative Equipment that allow users to work in a 100% Formaldehyde-Free Environment.

This is part of our range, and it is continuously growing with new high-quality products.



WARRANTY POLICY

- ✓ We guarantee, against any manufacturing defect, our supplied products **for 2 years**, from the as factory departure.
- ✓ This 2-year warranty applies as long as the equipment is used under normal conditions, in accordance with the applicable manual and the instructions for application, cleaning, repair, handling of possible defects and their continuous supervision.
- ✓ This warranty obligation consists **EXCLUSIVELY** in the replacement of the product, or part of it, recognized as defective, not including the costs of transport, installation, assembly or disassembly, or the damages that may have been caused.
- ✓ The product or component delivered in replacement of the defective one will be guaranteed in the same terms and conditions applicable to the replaced.
- ✓ The replaced product or component must be returned to **AnatHomic Solutions, S.L.U.**, carriage paid, within a maximum period of 15 days, from the shipment of the replaced products or components, after which **AnatHomic Solutions, S.L.U.** reserves the right to invoice replaced products or components.
- ✓ Products cannot be modified. They must be mounted as indicated in the assembly drawings. The settings cannot be altered.
- ✓ The anomalies produced by the transport must be informed "on site" and signed by the carrier. Otherwise, we cannot meet the warranty.
- ✓ Our entire range of manufactured products has a useful life of between 12 and 15 years, provided that the instructions for use and preventive maintenance described in the manuals of each equipment are followed.
- ✓ We guarantee the availability of spare parts: accessories, parts and pieces during the useful life of the equipment, provided that the instructions for use and preventive maintenance described in the manuals of each equipment are followed.
- ✓ Any warranty extension that may be required by the public or private sale process will be subject to negotiation and/or the presentation of specific documentation to be agreed between **AnatHomic Solutions, S.L.U.** and the applicant.

EXCLUSIONS

- ✓ The warranty will be void if full payment has not been made for the guaranteed products.
- ✓ The warranty does not cover damages caused by negligence in the transport, storage or handling of equipment, by operators or procedures external to **AnatHomic Solutions, S.L.U.**
- ✓ Likewise, the warranty will be voided if the maintenance, supervision, timely notification of possible defects, cleaning and repair operations specified in the user manual or instructions provided by **AnatHomic Solutions, S.L.U.** are not carried out.
- ✓ The warranty excludes breakdowns caused by incorrect or improper use of the equipment, or by defects caused by lack of maintenance or user training.
- ✓ In the same way, all responsibility and guarantee derived from the installation of the equipment is declined, if it has not been executed by authorized personnel or following the instructions indicated by **AnatHomic Solutions, S.L.U.**
- ✓ **AnatHomic Solutions, S.L.U.** declines all responsibility for delays justified by our suppliers such as: strikes, fires, catastrophes, etc., or modifications made once the order has been placed.
- ✓ The warranty and CE certification of our products is canceled, in the case of use of any type of spare part, part or consumable that is not provided and validated by **AnatHomic Solutions, S.L.U.**
- ✓ The warranty is void in the event of any type of alteration or technical modification or labeling of the equipment or its parts or consumables by any agent other than authorized by **AnatHomic Solutions, S.L.U.**



Scan our
vCard



C/del Ebre, 74 · 08223 Terrassa · Barcelona

T: +34 937 83 74 11 / +34 638 017 655

info@anathomicsolutions.com

www.anathomicsolutions.com

*Note: This document is for information purposes only. Any functionality, feature or specification are subject to change without notice. **AnatHomic Solutions SLU** reserves the right to improve this document or the product described without the obligation to communicate such changes.*

*Nota: Este documento tiene efectos meramente informativos. Cualquier funcionalidad, característica o especificación está sujeta a cambios sin previo aviso. **AnatHomic Solutions SLU** se reserva el derecho de mejorar este documento o el producto descrito, sin obligación de comunicar dichos cambios.*