

## MAXIMUM SECURITY



## FORMALDEHYDE DISPENSER

**Model DSF**

Cod. DSF-005-00-00



**AnatHomic Solutions**  
**PATHOLOGY LABORATORY**

## Why to choose our Formaldehyde Dispenser?

- ✓ Because it offers maximum safety and comfort to the user: it has a powerful rear and upper ventilation, and ventilation in the formaldehyde storage cubicle.
- ✓ Because it is possible to configure it according to the needs of users.
- ✓ Because its glossy finish facilitates optimal sanitization.
- ✓ Because the walls are covered with insulating material to reduce noise as much as possible.
- ✓ Because it has a filter counter that allows you to maximise your use.
- ✓ Because it is the only dispenser on the market that uses biodegradable filters that are very easy to replace.
- ✓ For its convenient replacement price of filters manufactured by AnatHomic Solutions.



## GENERAL CHARACTERISTICS


Stainless steel construction	AISI – 304
Sheet thicknesses	1.5 - 2 mm
Overall Dimensions	700 x 700 x (H)1975 mm
Connections: • Single-phase electrical connection 230V 50Hz + ground	
Electrical power	300 W
Control Panel	✓
Safety switch (current sectioning)	✓
Sides with window (safety glass)	3+3 laminated glass
Front window with safety glass	Laminated glass 4+4
Formaldehyde dispensing tap	✓
Vacuumed bottom housing for storage of formaldehyde tanks: • 1 tank of new 10L formaldehyde • 1 tank of used formaldehyde 25L • Tank filling and emptying sensors • Quick connect + anti-drip valves	✓
Automatic formaldehyde dispenser pedal	✓
Overhead LED Lighting	2
Operating Temperature	5°C – 45°C




## SUCTION SYSTEM

Double suction zone: front and top	✓
Centrifugal electric fan motors	1
Activated carbon filters impregnated with potassium permanganate	1
Filter life control	✓
First Notice Prior to Filter Change	800 h
End-of-life of filters	850 h
Suction speed adjustable from 10% to 100%	±5% intervals
Maximum suction capacity	1030 m3 / h

## SUCTION FLOW DIAGRAM

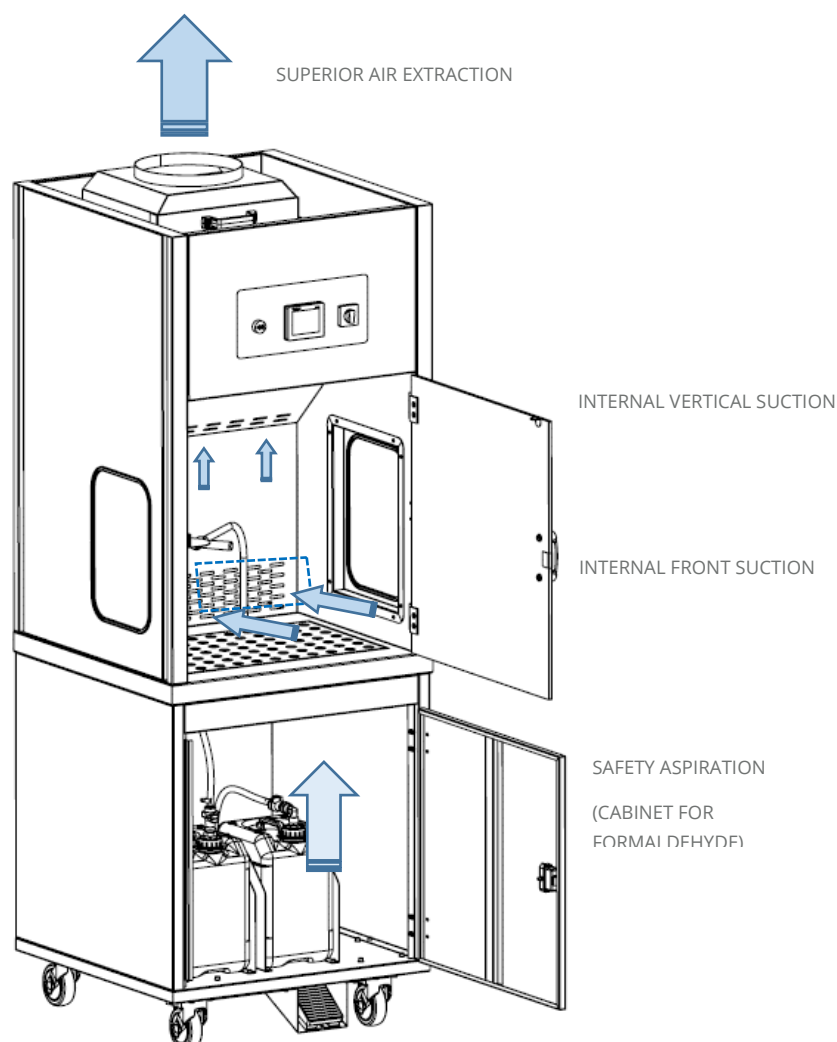


## DIMENSIONS

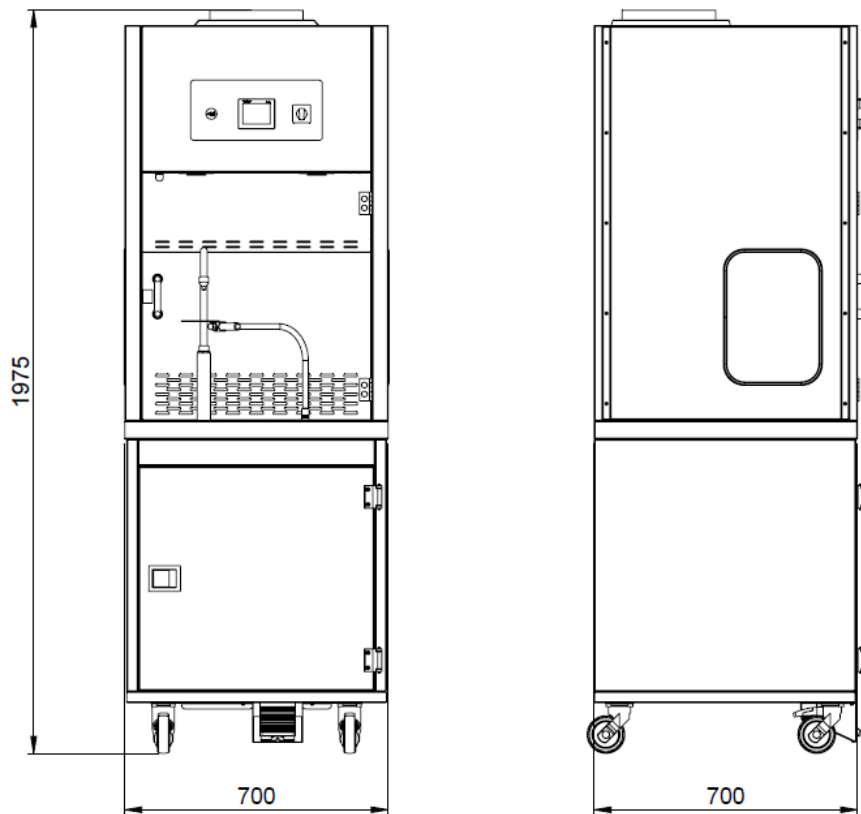

FORMALDEHYDE DISPENSER  
**DSF**

## OZONE DISINFECTION SYSTEM

All our stations come equipped with OZONIC WASH 4.0, which uses the power of ozone to eliminate pesticides, fungi, organic matter, pollutants, and viruses. In addition, it is 100% toxic-free, without adding pollutants to the water or the environment.

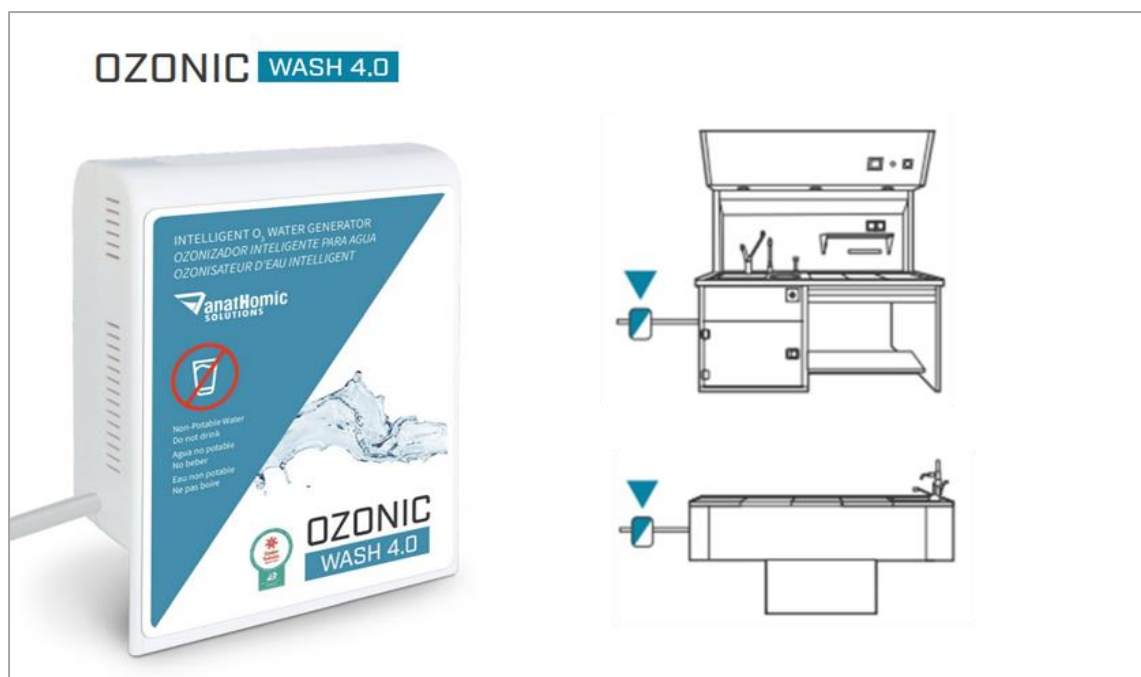

## BIODEGRADABLE FILTERS WITH POTASSIUM PERMANGANATE

All our stations are equipped with biodegradable filters made from corn and filled with potassium permanganate to disinfect exhaled particles into the environment. These filters are biodegradable and manufactured in-house, which allows us to offer the most competitive spare parts in terms of price and the generation of savings in the short term.



## PACKAGING

12-01 Phytosanitary Packaging + VCI Sleeve



## CERTIFICATIONS



UNE-EN ISO 9001:2015



UNE-EN ISO 13485:2018



UNE-EN ISO 14001:2015





Your Safety is our Motivation to Design and Manufacture Innovative Equipment that allows you to work in a 100% Formaldehyde-Free environment.

This is part of our continuously growing catalogue with new products that meet high quality standards.



## WARRANTY POLICY

- ✓ We guarantee, against any manufacturing defect, our supplied products **for 2 years**, from the date documented as leaving the factory.
- ✓ This 2-year warranty applies as long as the equipment is used under normal conditions, in accordance with the applicable manual and instructions for application, cleaning, repair, handling of possible defects and their continuous monitoring.
- ✓ This warranty obligation consists **EXCLUSIVELY** of the replacement of the product, or part thereof, recognized as defective, not including the costs of transport, installation, assembly or disassembly, or any damages that may have been caused.
- ✓ The product or component delivered in substitution of the defective one will be guaranteed under the same terms and conditions applicable to the substituted.
- ✓ The replaced product or component must be returned to **AnatHomic Solutions, S.L.U.**, postage paid, within a maximum period of 15 days, from the dispatch of the replaced products or components, after which **AnatHomic Solutions, S.L.U.** reserves the right to invoice replaced products or components.
- ✓ Products cannot be modified. They must be assembled as indicated in the assembly drawings. The settings cannot be altered.
- ✓ Anomalies caused by transport must be reported "on site" and signed by the carrier. Otherwise, we cannot honour the warranty.
- ✓ Our entire range of products have a useful life of between 12 and 15 years, as long as the instructions for use and preventive maintenance described in the manuals of each piece of equipment are followed.
- ✓ We guarantee the availability of spare parts: accessories, parts, and pieces for the useful life of the equipment, as long as the instructions for use and preventive maintenance described in the manuals of each equipment are followed.
- ✓ Any extension of the warranty that may be required by the public or private sale process will be subject to negotiation and/or the presentation of specific documentation to be agreed between **AnatHomic Solutions, S.L.U.** and the applicant.

## EXCLUSIONS

- ✓ The warranty will be void if full payment has not been made for the guaranteed products.
- ✓ The warranty does not cover damages caused by negligence in the transport, storage, or handling of the equipment, by operators or procedures external to **AnatHomic Solutions, S.L.U.**
- ✓ Likewise, the warranty will be void if the maintenance, supervision, timely notification of possible defects, cleaning and repair operations specified in the user manual or instructions provided by **AnatHomic Solutions, S.L.U.** are not carried out.
- ✓ The warranty excludes breakdowns caused by incorrect or improper use of the equipment, or by defects caused by lack of maintenance or user training.
- ✓ In the same way, all responsibility and warranty derived from the installation of the equipment is declined, if it has not been carried out by authorized personnel or following the instructions indicated by **AnatHomic Solutions, S.L.U.**
- ✓ **AnatHomic Solutions, S.L.U.** It declines all responsibility for delays justified by our suppliers such as: strikes, fires, disasters, etc., or modifications made after the order has been placed.
- ✓ The warranty and CE certification of our products is cancelled in the case of use of any type of spare part, part or consumable that is not provided and validated by **AnatHomic Solutions, S.L.U.**
- ✓ The warranty is void in the event of any type of alteration or technical modification or labelling of the equipment or its parts or consumables by any agent other than not authorised by **AnatHomic Solutions, S.L.U.**



Scan our  
vCard



C/del Ebre, 74 · 08223 Terrassa · BarcelonaT: +34 937 83 74 11 / +34 638 017 655  
info@anathomicsolutions.com

[www.anathomicsolutions.com](http://www.anathomicsolutions.com)

*Note: This document is for information purposes only. Any functionality, feature or specification are subject to change without notice. AnatHomic Solutions SLU reserves the right to improve this document, or the product described without the obligation to communicate such changes.*

*Note: This document is for informational purposes only. Any functionality, feature, or specification is subject to change without notice. AnatHomic Solutions SLU reserves the right to improve this document or the product described, without obligation to communicate such changes.*