

OUTSTANDING ERGONOMIC FOR TWO USERS



GROSSING STATION

Model ET-101-2

Code ET-XXX-XX-XX



AnatHomic Solutions
PATHOLOGY LABORATORY

Why to choose our ET-101?

- ✓ Because it offers great performance and ergonomics to the user, at a very competitive price using the highest quality in materials and in manufacturing process
- ✓ Because the level of security complies with the requirements to protect the user
- ✓ Because the walls are insulated to reduce the noise
- ✓ Because is the only one in the market that includes OZONIC wash system to disinfect water
- ✓ Because is the only one in the market that uses and includes biodegradable filters, very easy to be replaced
- ✓ Because of the very convenient pricing of filters manufactured by AnatHomic Solutions
- ✓ Because of easy maintenance, due to front and upper access
- ✓ Because the space on both sides of the Grossing Station is not used and instead other equipment, furniture or devices can be placed
- ✓ Because AnatHomic Solutions can provide, and assist in installation now or later, many complementary optional accessories with the same quality to make the most from our great Grossing Station: macrophotography system, dictaphone, shredder, water curtain, and many more.
- ✓ Because we are specialists in the manufacture and development of pathological anatomy equipment, and we know them so deeply that we are able to give the user exactly what he needs.
- ✓ Because we can offer very close and fast answer to our customer requirements, thanks to our great technical and commercial team.




GENERAL CHARACTERISTICS

Stainless steel construction	AISI – 304
Sheet metal thicknesses	1,5 - 2 mm
Overall dimensions	1500 x 950 x (h)1950 mm
Connections: <ul style="list-style-type: none"> • Single-phase electrical connection: 230V 50Hz + ground • Hot water supply 3/8" • Cold water supply 3/8" • Drain PVC Ø40 mm 	
Electrical power	300 W
Control panel <ul style="list-style-type: none"> • ON / OFF aspiration • ON / OFF lighting • Adjust ventilation speed percentage from 10% to 100% (intervals ±5 %) • Hour counter of filters • Visual and acoustic alarms with system notification (filters replacement, formaldehyde tanks levels) • Change language 	✓
Safety switch (power sectioning)	✓
Levelling feet	✓
Walls coated with insulating panels	✓
Formaldehyde dispenser tap	✓
Lower suction enclosure for storage formaldehyde tanks: <ul style="list-style-type: none"> • 1 new formaldehyde deposit 10L • 1 used formaldehyde deposit 10L • Filling and emptying level sensors • Quick connect valves + drip 	✓
Automatic formaldehyde dispenser pedal	✓
Washing sink	400 x 400 x 200 mm
Medical faucet (H/C water) by foot pedal	✓
Medical mixer shower with PVC cover <ul style="list-style-type: none"> • Integrated ozonated water disinfection system (H₂O+O₃) • Generation of ozone in water 0,2 – 1,0 ppm • Application: elimination of bacteria, viruses, mold, pesticides, and odours 	✓
Perforated removable trays	✓
Support shelf	✓
Waste bin with lid	✓
Magnetic instrument holder	✓
Auxiliary electrical outlets IP55	2
Top LED lighting	3
White LED lamp with flexible arm	15.000 lux a 50 cm
Foot support	✓

ASPIRATION SYSTEM





Aspiration from the lower, upper, and back area	✓
Centrifugal electro-fan motors	5
Activated carbon filters with permanganate impregnation	5
Filters useful life control system	✓
First notification for filters replacement	800 h
End of the filters service life	850 h
Adjustable suction speed 10% to 100%	Intervals of $\pm 5\%$
Maximum suction capacity	1404 m ³ / h

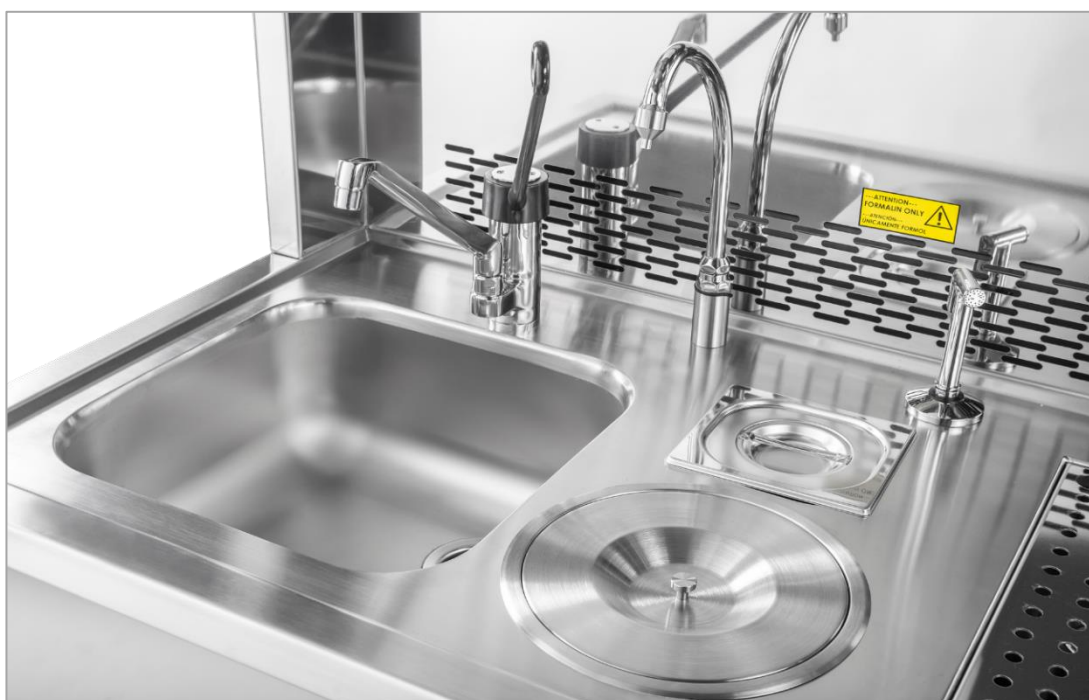

DIMENSIONS

ET-101-2

A	2400 mm	94.5 in
B	950 mm	37.4 in
C	1950 mm	76.8 in
D	850 mm	33.5 in
E	747 mm	29.4 in
F	420 mm	16.5
G	794 mm	31.3
H	728 mm	28.7

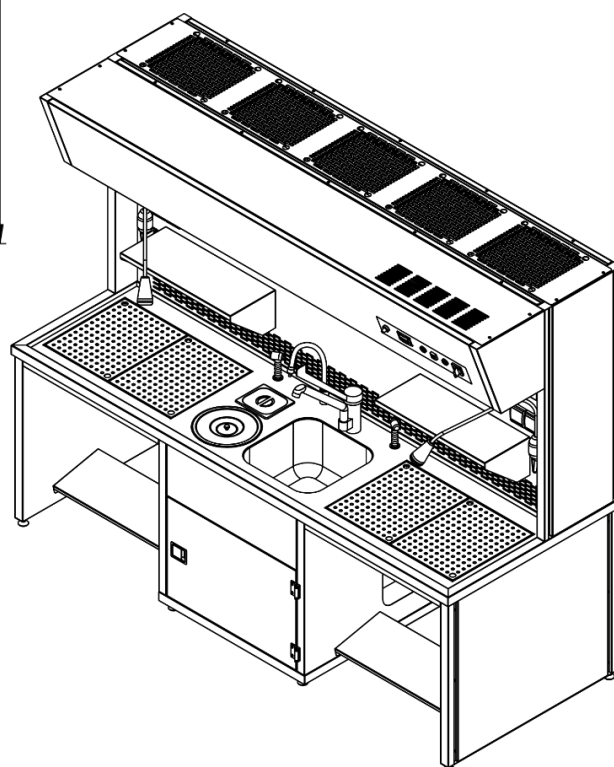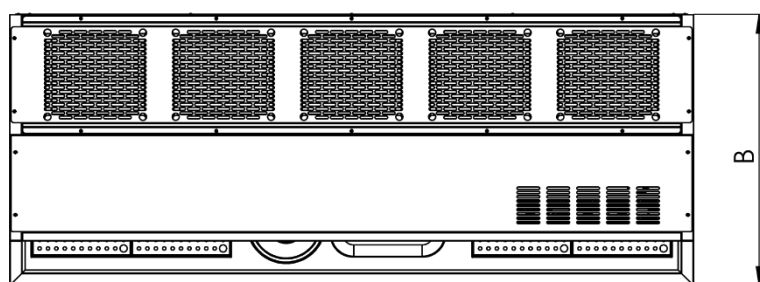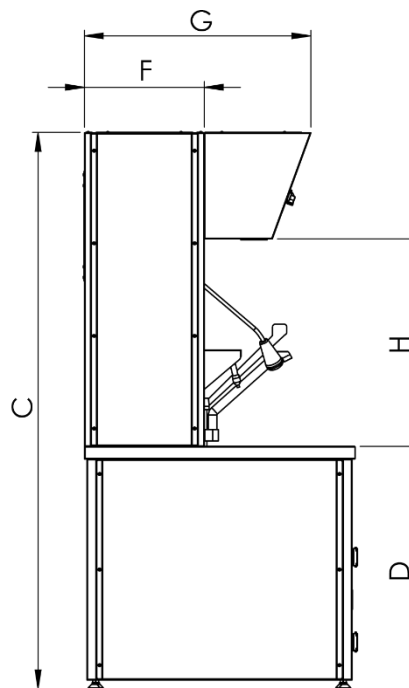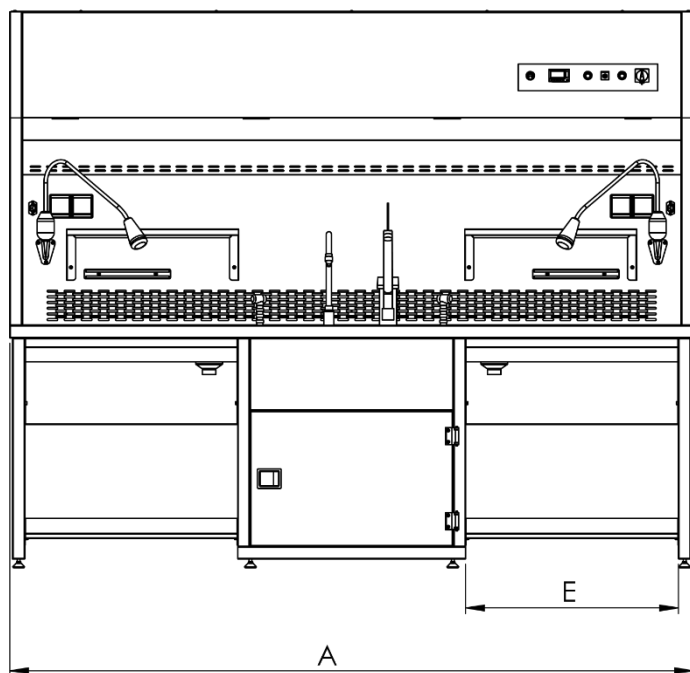

OZONE DISINFECTION SYSTEM

All our stations are equipped with OZONIC WASH 4.0, which through the power of ozone to eliminate pesticides, fungi, organic matter, pollutants, and viruses. In addition, it is 100% free of toxins, without adding contaminants to the water or the environment.



BIODEGRADABLE FILTERS WITH POTASSIUM PERMANGANATE


All our stations are equipped with biodegradable filters made from corn and filled with Potassium Permanganate to disinfect exhaled particles into the environment. These filters are manufactured by AnatHomic Solutions.



PACKAGING

12-01 Phytosanitary packaging + VCI cover

A	2550 mm	100.4 in
B	1100 mm	43.3 in
C	2200 mm	86.6 in
WHEIGHT	660 kg	1455 lb
Vol	6.17 m3	

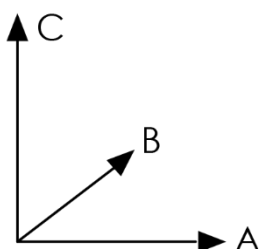

CERTIFICATIONS



UNE-EN ISO 9001:2015



UNE-EN ISO 13485:2018



UNE-EN ISO 14001:2015



Your Safety is our Motivation for Design and Manufacture Innovative Equipment that allow users to work in a 100% Formaldehyde-Free Environment.

This is part of our range, and it is continuously growing with new high-quality products.



WARRANTY POLICY

- ✓ We guarantee, against any manufacturing defect, our supplied products **for 2 years**, from the as factory departure.
- ✓ This 2-year warranty applies as long as the equipment is used under normal conditions, in accordance with the applicable manual and the instructions for application, cleaning, repair, handling of possible defects and their continuous supervision.
- ✓ This warranty obligation consists **EXCLUSIVELY** in the replacement of the product, or part of it, recognized as defective, not including the costs of transport, installation, assembly or disassembly, or the damages that may have been caused.
- ✓ The product or component delivered in replacement of the defective one will be guaranteed in the same terms and conditions applicable to the replaced.
- ✓ The replaced product or component must be returned to **AnatHomic Solutions, S.L.U.**, carriage paid, within a maximum period of 15 days, from the shipment of the replaced products or components, after which **AnatHomic Solutions, S.L.U.** reserves the right to invoice replaced products or components.
- ✓ Products cannot be modified. They must be mounted as indicated in the assembly drawings. The settings cannot be altered.
- ✓ The anomalies produced by the transport must be informed "on site" and signed by the carrier. Otherwise, we cannot meet the warranty.
- ✓ Our entire range of manufactured products has a useful life of between 12 and 15 years, provided that the instructions for use and preventive maintenance described in the manuals of each equipment are followed.
- ✓ We guarantee the availability of spare parts: accessories, parts and pieces during the useful life of the equipment, provided that the instructions for use and preventive maintenance described in the manuals of each equipment are followed.
- ✓ Any warranty extension that may be required by the public or private sale process will be subject to negotiation and/or the presentation of specific documentation to be agreed between **AnatHomic Solutions, S.L.U.** and the applicant.

EXCLUSIONS

- ✓ The warranty will be void if full payment has not been made for the guaranteed products.
- ✓ The warranty does not cover damages caused by negligence in the transport, storage or handling of equipment, by operators or procedures external to **AnatHomic Solutions, S.L.U.**
- ✓ Likewise, the warranty will be voided if the maintenance, supervision, timely notification of possible defects, cleaning and repair operations specified in the user manual or instructions provided by **AnatHomic Solutions, S.L.U.** are not carried out.
- ✓ The warranty excludes breakdowns caused by incorrect or improper use of the equipment, or by defects caused by lack of maintenance or user training.
- ✓ In the same way, all responsibility and guarantee derived from the installation of the equipment is declined, if it has not been executed by authorized personnel or following the instructions indicated by **AnatHomic Solutions, S.L.U.**
- ✓ **AnatHomic Solutions, S.L.U.** declines all responsibility for delays justified by our suppliers such as: strikes, fires, catastrophes, etc., or modifications made once the order has been placed.
- ✓ The warranty and CE certification of our products is canceled, in the case of use of any type of spare part, part or consumable that is not provided and validated by **AnatHomic Solutions, S.L.U.**
- ✓ The warranty is void in the event of any type of alteration or technical modification or labeling of the equipment or its parts or consumables by any agent other than authorized by **AnatHomic Solutions, S.L.U.**



Scan our
vCard



C/del Ebre, 74 · 08223 Terrassa · Barcelona

T: +34 937 83 74 11 / +34 638 017 655

info@anathomicsolutions.com

www.anathomicsolutions.com

*Note: This document is for information purposes only. Any functionality, feature or specification are subject to change without notice. **AnatHomic Solutions SLU** reserves the right to improve this document or the product described without the obligation to communicate such changes.*

*Nota: Este documento tiene efectos meramente informativos. Cualquier funcionalidad, característica o especificación está sujeta a cambios sin previo aviso. **AnatHomic Solutions SLU** se reserva el derecho de mejorar este documento o el producto descrito, sin obligación de comunicar dichos cambios.*